


Genotyping Report

Strain ID	T023018	Strain Type	CKO(Cas9)	Genetic Background	C57BL/6JGpt
Designer	Zifan Lin	Gene Name	Cops8		

1. Strategy of Genotyping

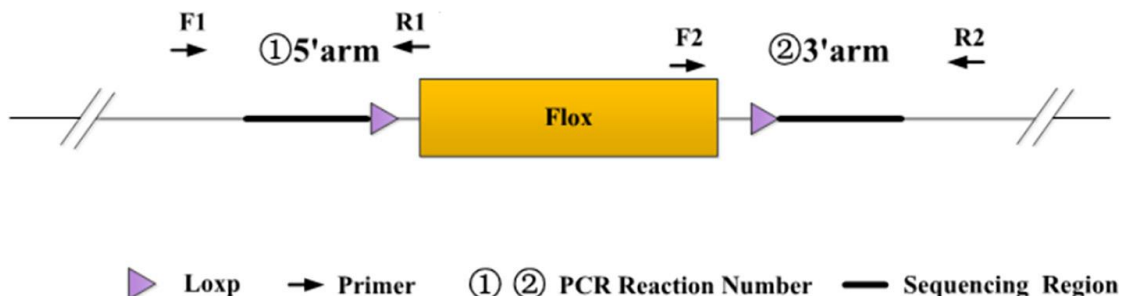

Wild type: ①PCR reaction obtains none band; ②PCR reaction obtains a single WT band.

Heterozygote: ①PCR reaction obtains a Targeted band; ②PCR reaction obtains a WT band and a Targeted band.

Homozygote: ①PCR reaction obtains a Targeted band; ②PCR reaction obtains a single Targeted band.

Note: The sizes of WT and Targeted band are shown below.

2. Primer Information

PCR No.	Primer No.	Sequence	Band Size
①(5'arm)	T023018-F1	CGGACAACGCCTTCAGTTTCAG	WT:0bp
	T023018-R1	CCAAGTACCTTGGGCAAGAACAT	Targeted:301bp
②(3'arm)	T023018-F2	GGCTTCTCTGCCATGATGCTGAT	WT:266bp
	T023018-R2	CCAGAAACATGAACAGCAGAGGTC	Targeted:372bp

3. Gel Image & Conclusion



Note: P:Positive control; WT: Wildtype control; B: Blank control (ddH₂O); M: DNA Ladder

4. PCR Condition

PCR Reaction Component			
Seg.	reaction component		Volume (μl)
1	2 × Taq Master Mix , Dye Plus, (Vazyme P112-03)		12.5
2	ddH2O		9.5
3	Primer A(10pmol/μl)		1
4	Primer B(10pmol/μl)		1
5	Template(≈100ng/μl)		1
PCR program			
Seg.	Temp.	Time	Cycle
1	95℃	5min	20×
2	98℃	30s	
3	65℃（-0.5℃/cycle）	30s	
4	72℃	45s	
5	98℃	30s	20×
6	55℃	30s	
7	72℃	45s	
8	72℃	5min	
9	10℃	hold	

Reviewer:Ting Sun

Date:2021-6-1