


Genotyping Report

Strain ID	T032377	Strain Type	KO(Cas9)	Genetic Background	C57BL/6JGpt
Designer	Zifan Lin	Gene Name	Cops7b		

1. Strategy of Genotyping

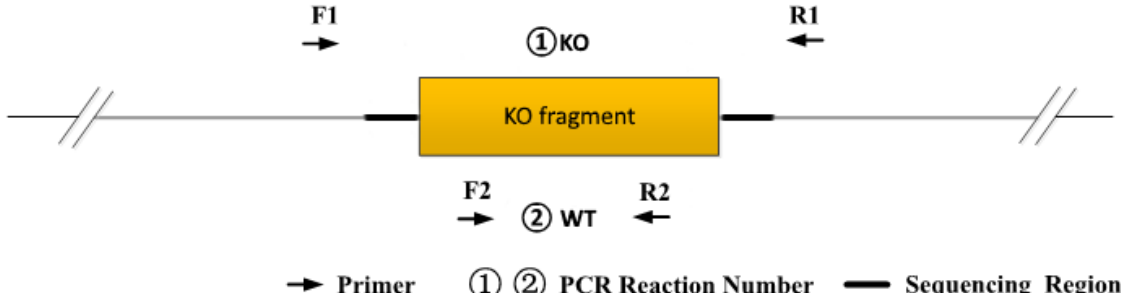


Wild type: ①PCR reaction obtains a single WT band; ②PCR reaction obtains a single WT band.
Heterozygote: ①PCR reaction obtains a WT band and a KO band; ②PCR reaction obtains a WT band.
Homozygote: ①PCR reaction obtains a single KO band; ②PCR reaction without product.
Note: 1)The sizes of WT and Targeted band are shown below.
2)If the WT band is too large, it may not be possible to obtain a WT band.

2. Primer Information

PCR No.	Primer No.	Sequence	Band Size
PCR①	T032377-F1	CAACCTCCACTTGTCACTTTAGATGC	WT:10201bp KO: 482bp
	T032377-R1	GCTCTGGTTTCATAATCTCCAAATCAC	
PCR②	T032377-F2	ATCCTTCCCACCTCCCTACACAC	WT: 427bp KO:0bp
	T032377-R2	GTGGTAGGTGACTCCAAGGAAAC	

3. Gel Image

cctcaggcttctgcctatgtaggatgg---9719bp---ttgtggggaagacaactgtgaataggcc

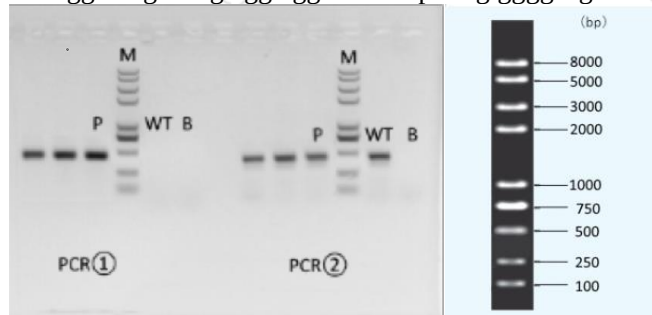

Note: P: Positive control; WT: Wildtype control; B: Blank control (ddH₂O); M: DNA Ladder

4. PCR Condition

PCR Reaction Component			
Seg.	reaction component		Volume (μl)
1	2 × Taq Master Mix , Dye Plus, (Vazyme P112-03)		12.5
2	ddH2O		9.5
3	Primer A(10pmol/μl)		1
4	Primer B(10pmol/μl)		1
5	Template(≈100ng/μl)		1
PCR program			
Seg.	Temp.	Time	Cycle
1	95℃	5min	20×
2	98℃	30s	
3	65℃（-0.5℃/cycle）	30s	
4	72℃	45s	
5	98℃	30s	20×
6	55℃	30s	
7	72℃	45s	
8	72℃	5min	
9	10℃	hold	

Reviewer:Ting Sun

Date:2021-6-16